CHANNEL CLEANING BRUSH

BW-18V

USA: CAUTION:
Federal law restricts this device to
sale by or on the order of a physician.

1

INTENDED USE

This instrument has been designed to be used for brushing the suction channels of OES TRANSDUODENAL CHOLEDOCHOFIBERSCOPE CHF-B20 and BRONCHOFIBERVIDEOSCOPE BF-MP190F/290F*1.
Do not use this product for any purpose other than its intended use.
*1 These products may not be available in some areas.

2

SYMBOLS

The meaning(s) of the symbol(s) shown on the back cover of this instruction manual are as follows:

Manufacturer

Authorized representative in the European Community

Importer (into European Union)

For US Customers only

For a Symbols Glossary, visit us: <http://www.olympus-global.com/en/common/pdf/symbolsglossary.pdf>

3

SIGNAL WORDS

The following signal words are used throughout this manual:

WARNING

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION


Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices or potential equipment damage.

NOTE

Indicates additional helpful information.

4

INSTRUMENT NOMENCLATURE

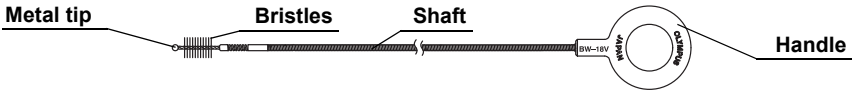5

PREPARATION, INSPECTION AND USE

CAUTION

- Before and after use, visually inspect the exterior of the Channel Cleaning Brush for breakage etc. In particular, confirm that Metal Tip is firmly attached. Using a damaged or worn brush may damage the endoscope being cleaned.
- The Channel Cleaning Brush is a consumable item. The bristle section can be damaged, by, for example, bending or buckling. If the bristle section becomes detached and remains inside the channel, insert (for example) biopsy forceps to retrieve it. Depending on the location, the bristle section may not be recoverable. In this case, contact Olympus. If the bristle section is left unrecovered within the scope, and the scope is used in this condition, the unrecovered portion may become dislodged from the scope during the patient procedure.
- To avoid splattering of detergent solution when the Channel Cleaning Brush is pulled out, keep the endoscope away from your face.

Before brushing, read the instruction manual for the compatible endoscope being cleaned thoroughly and perform the necessary preparation and inspection. Then, use the Channel Cleaning Brush according to the instruction manual for the compatible endoscope.

6

REPROCESSING

WARNING

- Patient debris and reprocessing chemicals are hazardous. Wear protective equipment.
- With the cleaning, disinfection and sterilization methods stated in this instruction manual, prions, which are considered to be the pathogenic substance of the Creutzfeldt-Jakob disease (CJD) cannot be destroyed or inactivated. When using this instrument on a patient with CJD or variant Creutzfeldt-Jakob disease (vCJD), be sure to use this product for such patient only and/or immediately dispose of this product after use in an appropriate manner. For methods to handle CJD, please follow the respective guidelines in your country.
- This instrument is not durable, or does not have sufficient durability against the respective methods stated in the guidelines of each country for destroying or inactivating prions. For information on the durability against each method, please contact Olympus. If cleaning, disinfection and sterilization methods not stated in this instruction manual are performed, Olympus cannot guarantee the effectiveness, safety and durability of this instrument. Make sure to confirm that there is no abnormality before use, and use under responsibility of a physician. Do not use if any abnormality is found.

CAUTION

When performing ultrasonic cleaning, disinfection or sterilization, do not coil the shaft with a diameter less than 15 cm. The shaft can be damaged if coiled too tightly.

6.1

Cleaning

- Soak the Channel Cleaning Brush in a low foaming detergent solution for the amount of time and at the temperature recommended by the detergent manufacturer.
- Clean the entire brush, while it is soaking in the low foaming detergent solution, using a sponge or lint-free cloth. In particular, clean off the bristle section with your fingertips to remove adhering debris.
- Rinse with water to remove the detergent solution completely.
- Perform ultrasonic cleaning, for 5 minutes at 38 – 47 KHz.
- Rinse in water, then thoroughly dry with lint-free cloth.

CAUTION

If the brush is not thoroughly cleaned, residual organic debris may solidify and it may be difficult to effectively disinfect or sterilize the brush.

6.2

Disinfection

- After cleaning the brush in accordance with Section 6.1, “Cleaning” prepare a 2.0 – 3.2% glutaraldehyde solution. Follow the manufacturer’s instructions regarding its concentration and temperature.
- Soak the Channel Cleaning Brush in a basin filled with disinfectant solution, then remove bubbles from the outer surface, using a lint-free cloth. In particular, rub the bristle section with your fingers to remove bubbles.
- Soak in disinfectant solution for the amount of time recommended by the disinfectant manufacturer.
- Remove from the disinfectant solution, then rinse thoroughly with sterile water.
- Thoroughly dry with a sterile lint-free cloth.

NOTE

Contact Olympus for the names of specific disinfectants with which the Channel Cleaning Brush is compatible.

6.3

Sterilization

- After cleaning the brush in accordance with Section 6.1, “Cleaning”, seal the brush in packaging appropriate for steam (autoclave) sterilization according to your hospital’s protocol.
- Autoclave the brush (steam sterilization) under the following conditions. Operate the autoclave in accordance with its instruction manual.
- Utilize the autoclave drying cycle, or leave the brush in the autoclave at the conclusion of the cycle with the door ajar, to ensure that the sterile packs are dry before handling.

Sterilizer	Temperature	Exposure time
Prevacuum	132 – 134°C (270 – 274°F)	5 minutes

CAUTION

When autoclaving, make sure that the bristles are not crushed flat.

7

STORAGE

Store the Channel Cleaning Brush at a room temperature, in a dry, clean location away from direct sunlight.

CAUTION

- Do not coil the shaft with a diameter less than 15 cm. The shaft can be damaged if coiled too tightly.
- Do not store the Channel Cleaning Brush in the carrying case. Use the carrying case only for shipping the endoscope. Routinely storing the Channel Cleaning Brush in a humid, dark, non-ventilated environment such as the carrying case can present an infection control risk.

If you have any questions or comments about any information in this manual, contact Olympus.

Manufactured by

OLYMPUS MEDICAL SYSTEMS CORP.
2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan
Fax: (042)646-2429 Telephone: (042)642-2111

OLYMPUS AMERICA INC.
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.
Fax: (484)896-7128 Telephone: (484)896-5000

OLYMPUS KOREA CO., LTD.
(Seocho-dong, Majestat City Tower One), 2F & 3F, 12, Seocho-daero 38-gil, Seocho-gu, Seoul, 06655, Republic of Korea
Fax: (02)6255-3210 Telephone: 1544-3200

OLYMPUS SINGAPORE PTE LTD
438B Alexandra Road #03-07/12, Alexandra Technopark Bld B, Singapore 119968
Fax: 6834-2438 Telephone: 6834-0010

Distributed by

OLYMPUS EUROPA SE & CO. KG
(Premises/Goods delivery) Wendenstrasse 20, 20097 Hamburg, Germany
(Letters) Postfach 10 49 08, 20034 Hamburg, Germany
Fax: (040)23773-4656 Telephone: (040)23773-0

OLYMPUS LATIN AMERICA, INC.
5301 Blue Lagoon Drive, Suite 290 Miami, FL 33126-2097, U.S.A.
Fax: (305)261-4421 Telephone: (305)266-2332

OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.
A8F, Ping An International Financial Center, No. 1-3, Xinyuan South Road, Chaoyang District, Beijing, 100027 P.R.C.
Fax: (86)10-5976-1299 Telephone: (86)10-5819-9000

OLYMPUS AUSTRALIA PTY LTD
3 Acacia Place, Notting Hill, VIC 3168, Australia
Fax: (03)9543-1350 Telephone: (03)9265-5400

KEYMED (MEDICAL & INDUSTRIAL EQUIPMENT) LTD.
KeyMed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, United Kingdom
Fax: (01702)465677 Telephone: (01702)616333

OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
Elektrozavodskaya str. 27 bld 8, 107023 Moscow, Russia
Fax: (7)495-926-7072 Telephone: (7)495-926-7077

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